


BCA Assay w/ Cytation 5 (562 nm)

Source: doc_38f142c704f3_20260915-ES421M-ES477M-ES495M-colon-back-wing.xlsx | Date: 2026-09-15 15:22:52 | Reader: Cytation5 (S/N: 1710247)
Plate: 96 WELL PLATE | Temperature: 24.8 °C | Wavelength: 562 nm (Absorbance endpoint)

1. 4PL Calibration Curve & Parameters

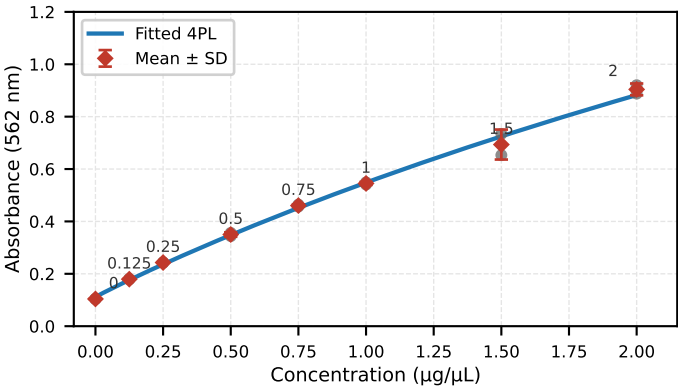

4PL Model Equation:
$$y = \text{bottom} + \frac{\text{top} - \text{bottom}}{1 + (\text{EC}_{50}/x)^{\text{HillSlope}}}$$


Fitted Parameters:

- Bottom = 0.111444
- Top = 4.0000
- EC50 = 8.49012 µg/µL
- Hill Slope = 0.96631

Goodness of Fit:


- R² = 0.996812
- RMSE = 0.014373 OD
- Standards: Wells A1:H2 duplicate series (2.0 to 0.0 µg/µL)
- Monotonic: Yes | Valid fit: Yes

2. Raw Absorbance Plate (OD 562 nm)

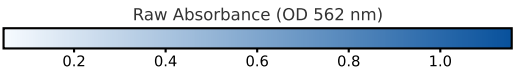

	1	2	3	4	5	6	7	8	9	10	11	12
A	0.888	0.92	0.865	0.937	0.674	1.156	0.637	0.553	0.539	0.595	0.484	0.05
B	0.734	0.653	0.439	0.475	0.324	0.516	0.352	0.281	0.285	0.313	0.271	0.049
C	0.541	0.549	0.048	0.048	0.048	0.048	0.047	0.048	0.05	0.049	0.05	0.047
D	0.457	0.464	0.048	0.048	0.047	0.048	0.049	0.049	0.048	0.055	0.048	0.048
E	0.344	0.356	0.048	0.049	0.048	0.048	0.046	0.048	0.048	0.047	0.048	0.048
F	0.242	0.244	0.05	0.048	0.049	0.048	0.047	0.049	0.048	0.05	0.048	0.048
G	0.179	0.18	0.048	0.049	0.048	0.047	0.049	0.047	0.048	0.049	0.048	0.047
H	0.104	0.104	0.048	0.05	0.048	0.048	0.051	0.05	0.049	0.049	0.048	0.048

3. Estimated Concentrations (µg/µL, Inverse 4PL)

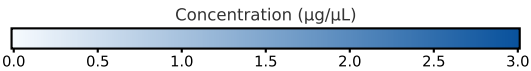

	1	2	3	4	5	6	7	8	9	10	11	12
A	2.019	2.127	1.942	2.186 △	1.350	3.011 △	1.244	1.012	0.975	1.126	0.832	
B	1.528	1.289	0.719	0.809	0.444	0.915	0.509	0.348	0.356	0.419	0.325	
C	0.980	1.002										
D	0.764	0.781										
E	0.491	0.519										
F	0.262	0.267										
G	0.130	0.132										
H	-0.013	-0.013										

Notes & Quality Control:

- Standards: Wells A1:H2 duplicate series (2.000 to 0.000 µg/µL BSA, dotted black border).
- △ Flagged Wells: Wells marked with △ (red borders) exceed the 2.0 µg/µL calibration range or detector OVRFLW (>4.0 OD).
- Diagnostics: Negative values are signed fit diagnostics, not physical concentrations. Blank wells indicate baseline (<0.065 OD).
- Recommendation: Samples with △ (e.g. A4, A6) should be diluted (e.g. 1:5 or 1:10) and re-assayed.

